



Filter Pipette Tips

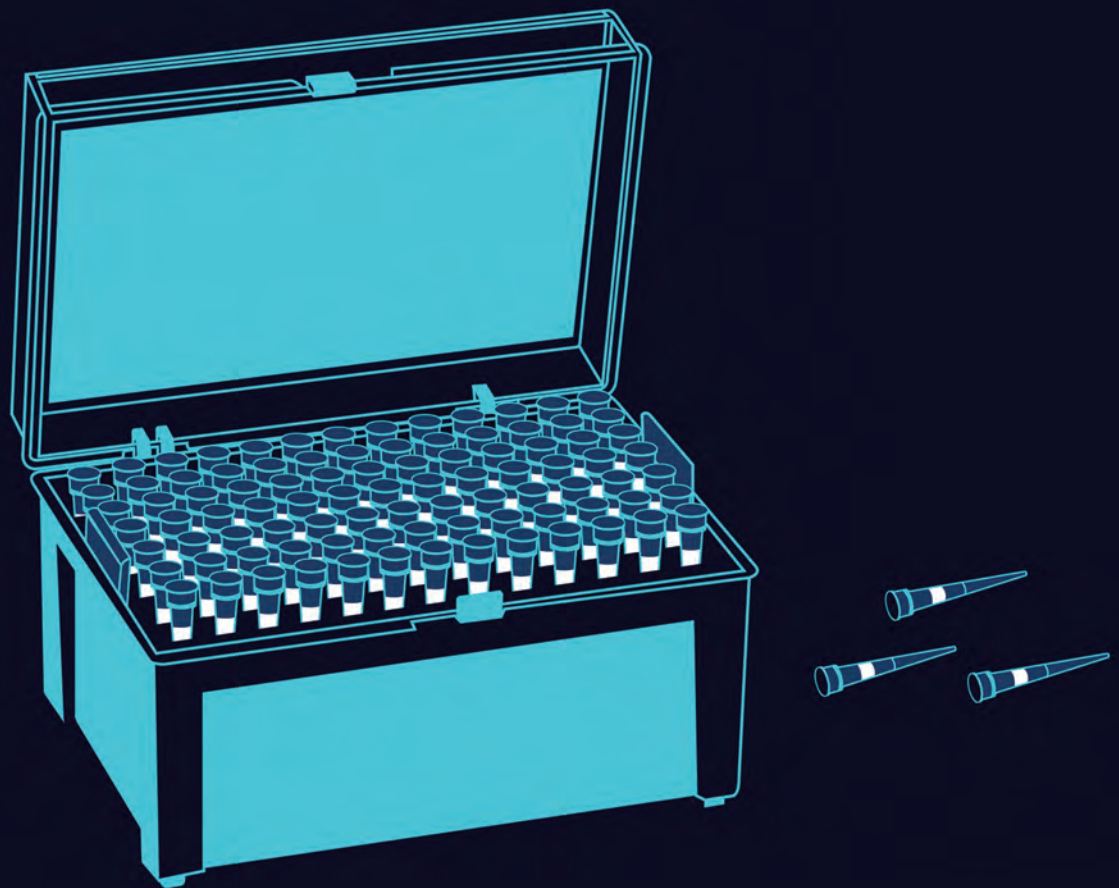

Your Trusted Partner in Life Sciences

www.abdoslabware.com

- Abdos filter pipette tips are manufactured from ultra pure USP class VI certified virgin polypropylene to ensure the highest quality & conforming to US FDA 21 CFR. Free from heavy metals, natural rubber & inhibitory plasticizers.
- Proprietary polyethylene super filter provides aerosol barrier and aqueous liquid by-pass prevention in over-pipetting hereby preventing cross contamination.
- Best suited for biological application.

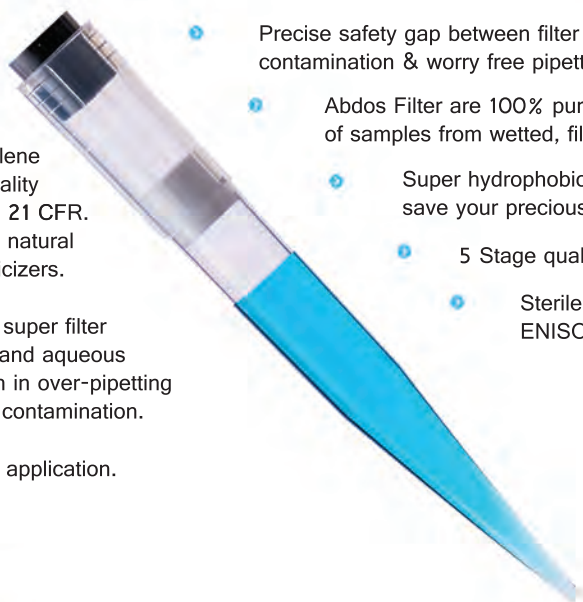

Precise safety gap between filter and liquid, filter and pipette shaft to avoid contamination & worry free pipetting.

- Abdos Filter are 100% pure hydrophobic polyethylene (PE) makes recovery of samples from wetted, filtered without contamination.
- Super hydrophobic filter that does not absorb samples, thus save your precious sample loss in pipetting error.
- 5 Stage quality control procedures ensure 100% quality & performance.
- Sterile tips has been validated according to ENISO 11137-2; 2012 for SAL 10^{-6} .



Filters used in Abdos filter pipette tips are non cytotoxic (ISO 10993-5 MRM, minimal essential media elution, GLP compliant). Abdos filters provide 99.9% bacterial aerosol filtration efficiency (BFE), (ASTM F21012), are non hemolytic (modified ASTM F 756-08 hemolysis direct contact GLP with human blood), inert and hydrophobic filter does not trap or block the liquid.



Code

Product Description

No. Per Case

P11021	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 0.2 - 10 μ l	10 Racks x 96 pcs (960 pcs)
P11023	ABDOS Filter Pipette Tips, Racked Sterile, Extended Graduated, 0.2 - 10 μ l	10 Racks x 96 pcs (960 pcs)
P11025	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 2 - 20 μ l	10 Racks x 96 pcs (960 pcs)
P11027	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 2 - 50 μ l	10 Racks x 96 pcs (960 pcs)
P11029	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 2 - 100 μ l	10 Racks x 96 pcs (960 pcs)
P11033	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 2 - 200 μ l	10 Racks x 96 pcs (960 pcs)
P11035	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 2 - 300 μ l	10 Racks x 96 pcs (960 pcs)
P11037	ABDOS Filter Pipette Tips, Racked Sterile, 100 - 1000 μ l	10 Racks x 96 pcs (960 pcs)

* For additional product or technical information, please visit www.abdoslabware.com. For distributor information contact labtech@abdosindia.com.

Abdos Labtech Pvt Ltd

720 - 722, DLF Tower A,
Jasola, New Delhi - 110025, India.

T: +91-11-44800980-986

F: +91-11-44800987

E: labtech@abdosindia.com

Europe Office

T: +44 - 7958295643

E: eu.labtech@abdosindia.com