

# Technical Data Sheet

Doc No – Ab/TDS/RB-44-II

Rev No : 01

## TECHNICAL DATA SHEET

**Abdos Labtech Pvt. Ltd.**

**Jumbo Test Tube Rack, RPP**

**Valid for Cat. No. P20443, P20444,**

### **1. Description/Specification:**

#### **1.1 Description:**

Jumbo Test Tube Rack, RPP

1.1.1 These innovative jumbo test tube racks are used for high capacity storage for tubes, vessels & vials.

#### **1.2 Design:**

Ergonomic design for easy handling Holes are keyed with molded-in- numbers and letters for easy identification.

#### **1.3 Performance Specifications:**

<b>Code</b>	<b>Tubes (mm)</b>	<b>Places</b>	<b>L x W x H (mm)</b>	<b>No. Per Case</b>
<b>P20443</b>	10-13	216	337x267x91	1
<b>P20444</b>	13-16	120	337x267x91	1

#### **1.4 Material:**

1.4.1 Complies with requirement of European Pharmacopoeia, 7<sup>th</sup> Edition (2011) & its supplements 7.5 (07/2012) monograph 3.2.2 Plastic container and closures for Pharmaceutical use.

1.4.2 The product has Passed the USP testing including class VI tests & has been assigned the FDA Drug master file DMF 12123

1.4.3 Material has Successful Passed the biological tests according to ISO 10993 external Communication devices for indirect blood contact for Prolonged period.

**1.5 Color:** White, Blue



Shiv Ganga Industrial Estate,  
Plot no. 1 Lakeshwari, Bhagwanpur,  
Roorkee - 247667 Uttarakhand, India  
**An ISO 9001:2008 Company**  
A Division of Abdos

**Regd. Office :**  
13A, Government Place (East)  
Kolkata – 700069, India  
**CIN :** U85100WB2008PTC123903

**P :** +91 (0) 1332 235277  
**E :** info@abdoslabware.com  
**W :** www.abdoslabware.com

# Technical Data Sheet

## 1.6 Sterilization:

**Sterile catalog number** from this lot have been pre-sterilized and dosimetrically released based on ANSI/AAMI/ISO 11137 *Sterilization of healthcare products-Requirements for validation and routine control- Radiation sterilization*. The product has been sterilized with Co-60 Gamma Radiation.

## 1.7 Quality control:

DNase, RNase, and Pyrogen Free

Testing of various characteristics in accordance with the valid Specification

## 1.8 Other information

Reusable, can be used in wide temperature range

## 1.9. Basic Features


**1.9.1 Temperature range:** -20<sup>0</sup>C to +121<sup>0</sup>C


**1.9.2 Autoclavability:** Yes, Autoclavable at 15 psi and 121<sup>0</sup>C for 20 minutes

**1.9.3 Chemical resistance:** See website.

## Autoclaving Guideline for Plastics

All plastic consumables are autoclavable at 121<sup>0</sup>C for 20 minutes, 15 psi/1atm. Do not use a dry cycle as this may cause deformation or weakening of the plastic. Allow temperature within the autoclave to return to at least 50<sup>0</sup>C before removing product.

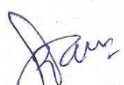  
Quality Assurance Incharge  
Abdos Labtech Pvt. Ltd.

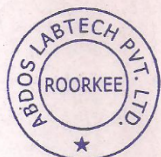

Shiv Ganga Industrial Estate,  
Plot no. 1 Lakeshwari, Bhagwanpur,  
Roorkee - 247667 Uttarakhand, India  
**An ISO 9001:2008 Company**  
A Division of Abdos

**Regd. Office :**  
13A, Government Place (East)  
Kolkata – 700069, India  
**CIN :** U85100WB2008PTC123903

**P :** +91 (0) 1332 235277  
**E :** info@abdoslabware.com  
**W :** www.abdoslabware.com