


ABDOS
LIFE SCIENCE

10 YEARS **ABDOS**
LIFE SCIENCE

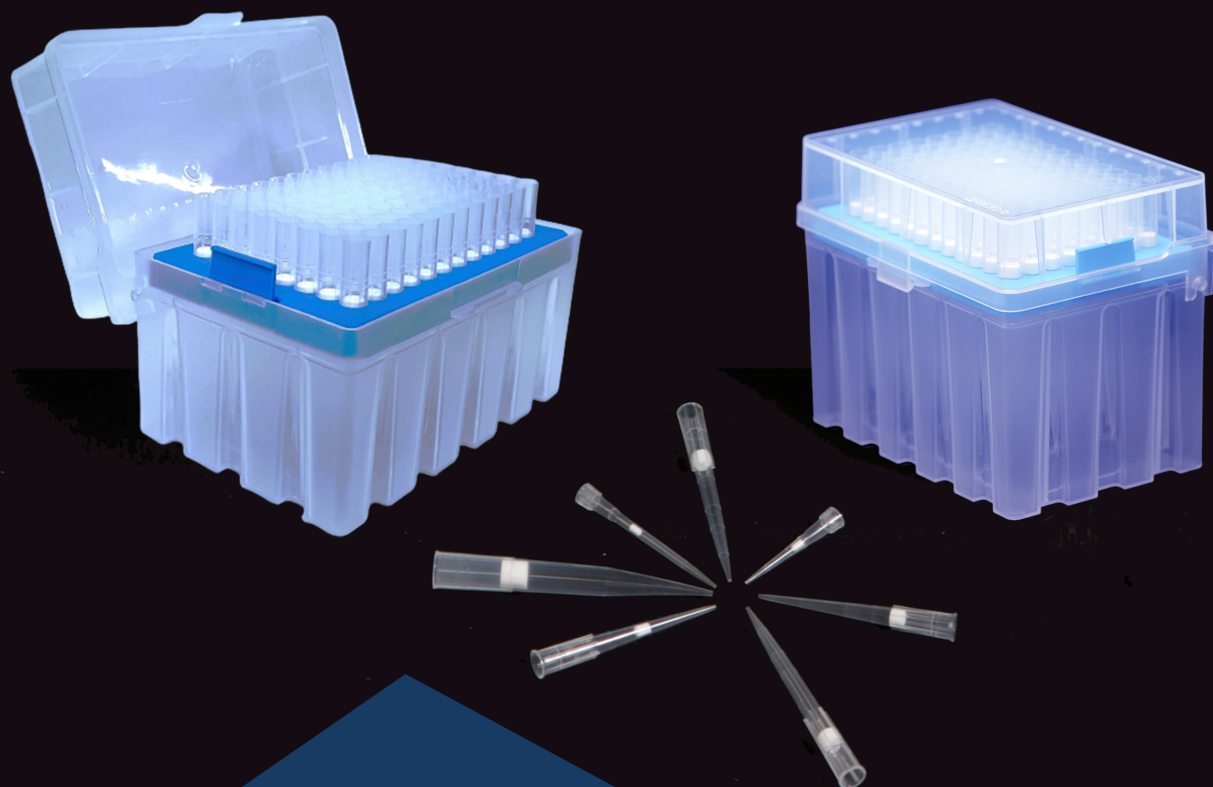

Filter Pipette Tips


Creating For Life Sciences

www.abdoslifesciences.com

Abdos Super Filter Prevents Both Liquid & Aerosol Passing Through The Filter.

Competitor Filter pipette Tips

Abdos Super Filter Pipette Tips

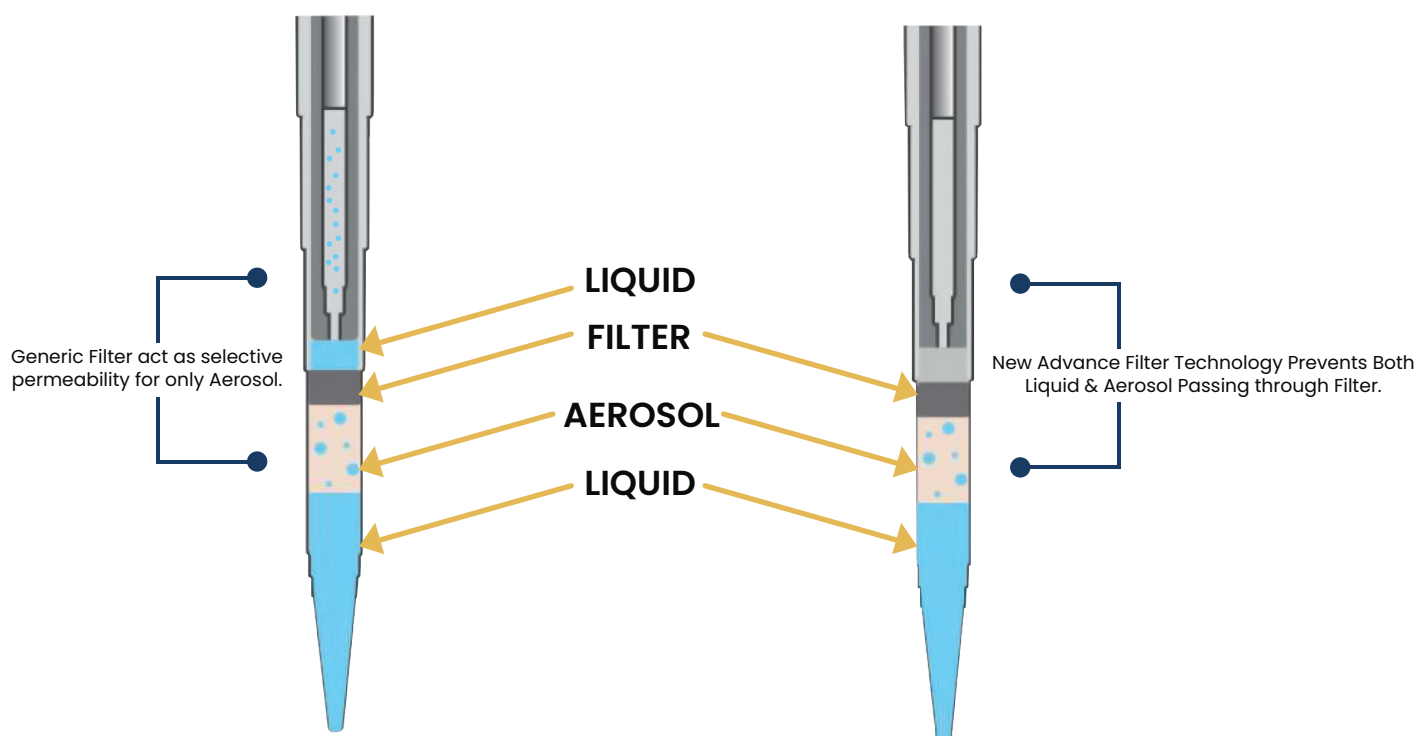


Abdos Super Filter Technology provides aerosol barrier & aqueous liquid* by-pass prevention in over-pipetting thus safeguarding precious samples against contamination.

Ensure Smooth & Hassle Free Pipetting with Increased Air Gap to Prevent Cross Contamination.

Abdos 200 μ l Filter Tips

INCREASED AIR GAP


Competitor 200 μ l Filter Tips

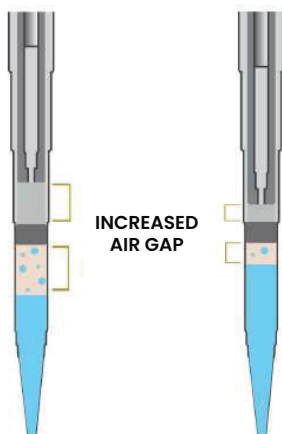

Experience The Finest Super Filter Technology With Abdos Filter Tips



- Universal fit, leak proof fitment.
- Precise safety gap between filter and liquid, filter and pipette shaft to avoid contamination & worry free pipetting.
- Abdos Filter is 100% pure porous polyethylene (PE) helps us to recover samples from wetted filter and does not trap or block the liquid into the filter.
- 5 Stage quality control procedures ensures 100% quality & performance.
- Inert hydrophobic filter does not trap or block the liquid, thus save your precious sample loss in pipetting error.
- Sterile tips has been validated according to ENISO 11137-2; 2012 for SAL 10^{-6} .
- Abdos filter pipette tips are manufactured from ultra pure USP class VI certified virgin polypropylene to ensure the highest quality & confirming to US FDA 21 CFR. Free from heavy metals, natural rubber & inhibitory plasticizers.
- Proprietary polyethylene super filter provides aerosol barrier and aqueous liquid by-pass prevention in over-pipetting hereby preventing cross contamination.
- Best suited for biological application.

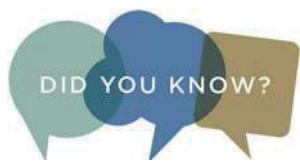

Filters used in Abdos Last Drop filter pipette tips are non cytotoxic (ISO 10993-5 MRM, minimal essential media elution, GLP compliant). Abdos filters provide 99.9% bacterial aerosol filtration efficiency (BFE),(ASTM F21012), are non hemolytic (modified ASTM F 756-08 hemolytic direct contact GLP with human blood), inert and hydrophobic filter that does not trap or block the liquid.

Code	Product Description	No.Per Case
P11021	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 0.1 - 10 µl	10 Racks x 96 pcs (960 pcs)
P11023	ABDOS Filter Pipette Tips, Racked Sterile, Extended Graduated, 0.1 - 10 µl	10 Racks x 96 pcs (960 pcs)
P11025	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 1 - 20 µl	10 Racks x 96 pcs (960 pcs)
P11027	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 1 - 50 µl	10 Racks x 96 pcs (960 pcs)
P11029	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 1 - 100 µl	10 Racks x 96 pcs (960 pcs)
P11033	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 1 - 200 µl	10 Racks x 96 pcs (960 pcs)
P11035	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 1 - 300 µl	10 Racks x 96 pcs (960 pcs)
P11037	ABDOS Filter Pipette Tips, Racked Sterile, 100 - 1000 µl	10 Racks x 96 pcs (960 pcs)
P11047	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 100 - 1000 µl	10 Racks x 96 pcs (960 pcs)
P11048	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 100 - 1250 µl	10 Racks x 96 pcs (960 pcs)
P10172	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 100 - 5000 µl	6 Racks x 24 pcs
P10173	ABDOS Filter Pipette Tips, Racked Sterile, Graduated, 500 - 10000 µl	6 Racks x 24 pcs
P11020	ABDOS Filter Pipette Tips, Racked, Graduated, 0.1 - 10 µl	10 Racks x 96 pcs (960 pcs)
P11022	ABDOS Filter Pipette Tips, Racked Extended Graduated, 0.1 - 10 µl	10 Racks x 96 pcs (960 pcs)
P11024	ABDOS Filter Pipette Tips, Racked, Bevelled Graduated, 1 - 20 µl	10 Racks x 96 pcs (960 pcs)
P11026	ABDOS Filter Pipette Tips, Racked, Bevelled Graduated, 1 - 50 µl	10 Racks x 96 pcs (960 pcs)
P11028	ABDOS Filter Pipette Tips, Racked, Bevelled Graduated, 1 - 100 µl	10 Racks x 96 pcs (960 pcs)
P11032	ABDOS Filter Pipette Tips, Racked, Graduated, 1 - 200 µl	10 Racks x 96 pcs (960 pcs)
P11034	ABDOS Filter Pipette Tips, Racked, 100 - 300 µl	10 Racks x 96 pcs (960 pcs)
P11036	ABDOS Filter Pipette Tips, Racked, 100 - 1000 µl	10 Racks x 96 pcs (960 pcs)
P11038	ABDOS Filter Pipette Tips, Racked, Graduated, 100 - 1000 µl	10 Racks x 96 pcs (960 pcs)
P11039	ABDOS Filter Pipette Tips, Racked, Graduated, 100 - 1250 µl	10 Racks x 96 pcs (960 pcs)
P11030	ABDOS Filter Pipette Tips, Racked, Graduated, 100 - 5000 µl	6 Racks x 24 pcs
P11031	ABDOS Filter Pipette Tips, Racked, Graduated, 500 - 10000 µl	6 Racks x 24 pcs




Filter Pipette Tips, Bulk in Resealable Bags



Filter Pipette Tips, Racked



Filter Pipette Tips, Sterile Racked

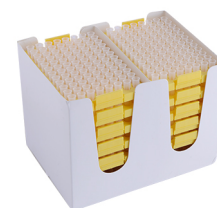

Filter Pipette Tips, Refill Pack