

# Technical Data Sheet

Doc No – Ab/TDS/PCR-104

Rev No: 07 Date 04/2023-24

## TECHNICAL DATA SHEET

**Abdos Labtech Pvt. Ltd.**

**AmpQlear™ PCR Strip Tubes without Cap 8 Strips, PP**

**Valid for Cat. No. P10706, P10708, P10732, P10733**

**SRN (Single Registration Number): IN-MF-000022272**

### 1. Description/Specification

#### 1.1 Description:

Strip Tubes without Cap 8 Strips, PP, DNase, RNase, Human DNA, ATP, Pyrogen Free, & PCR Inhibitors Free, IVD, CE & ISO18385 certified.

#### 1.2 Design:

Strip Tubes without Cap & with Cap

#### 1.3 Performance Specification:

Code	Capacity (ml)	Type	Height (mm)	No. per Case
<b>P10706</b>	0.2	8 Strip Standard Tubes	20.70	125 Strips
<b>P10708</b>	0.1	8 Strip Low profile Tubes	16.10	125 Strips
<b>P10732</b>	0.2	White 8 Strip Standard Tubes	20.70	125 Strips
<b>P10733</b>	0.1	White 8 Strip Low Profile Tubes	16.10	125 Strips

#### 1.4 Material:

1.4.1 Polypropylene (PP), Hydrophobic, confirm US FDA 21CFR

1.4.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for Pharmaceutical use.

1.4.3 The product has passed the USP testing including class VI tests ,Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.4.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE ,2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

1.4.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

## 1.5 Color:

Translucent (Natural)/White

## 1.6 Quality control:

Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release

## 1.7 Other information

1.7.1 DNase, RNase, DNA, Endotoxin and PCR inhibitor free.

1.7.2 Recommended for Single use only

1.7.3 Thin and uniform wall thickness to ensure optimal heat transfer reducing insulating effect.

1.7.4 Snap cap ensures tight fit eliminating risk of sample evaporation. Flat caps are optically clear for real time qPCR.

## 1.8 Basic Features

**1.8.1 Temperature range:** -80 °C to + 121°C

**1.8.2 Autoclavability:** Yes, Autoclavable at 15 psi and 121 °C for 20 minutes

**1.8.3 Manufactured in Aseptic Environment**

**1.8.4 Chemical resistance:** see website ([www.abdoslifesciences.com](http://www.abdoslifesciences.com))

**1.8.5 Tubes have ultra thin walls to reduce the insulating effect on the strips**

**1.8.6 Strips are available in both Standard & Low Profile wells**

## Autoclaving Guideline for Plastics

All plastic consumable are Autoclavable at 121°C for 20 minutes, 15 psi/1atm. Do not use dry cycle as this may cause deformation or weakening of plastic. Allow temperature within the autoclave to return to at least 50°C before removing product.

## Instruction for Use

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email: [labtech@abdosindia.com](mailto:labtech@abdosindia.com)  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet


## WARNING

**DO NOT** store strong oxidizing agents.

**DO NOT** place any plastic labware in a flame.

**DO NOT** mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids.

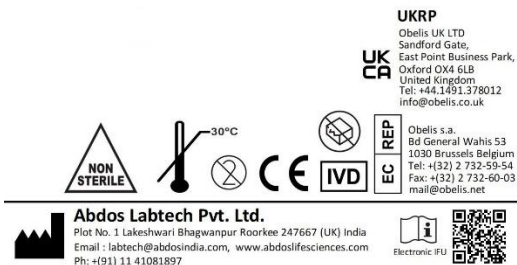

Quality Manager

Abdos Labtech Private Limited

Employee Code 30110



**Abdos Labtech Private Limited**  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheshwari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903

**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

**Ph. No: +91- (0) 1332235277**  
**Email:labtech@abdosindia.com**  
**Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)**  
**GST: 05AAGCA8605A1Z7**